


For Small & Medium sized hospital and Nursing hospital



RAK-U06M2E

/RAK-U03M2E

CE 0068
CE MDD Certified

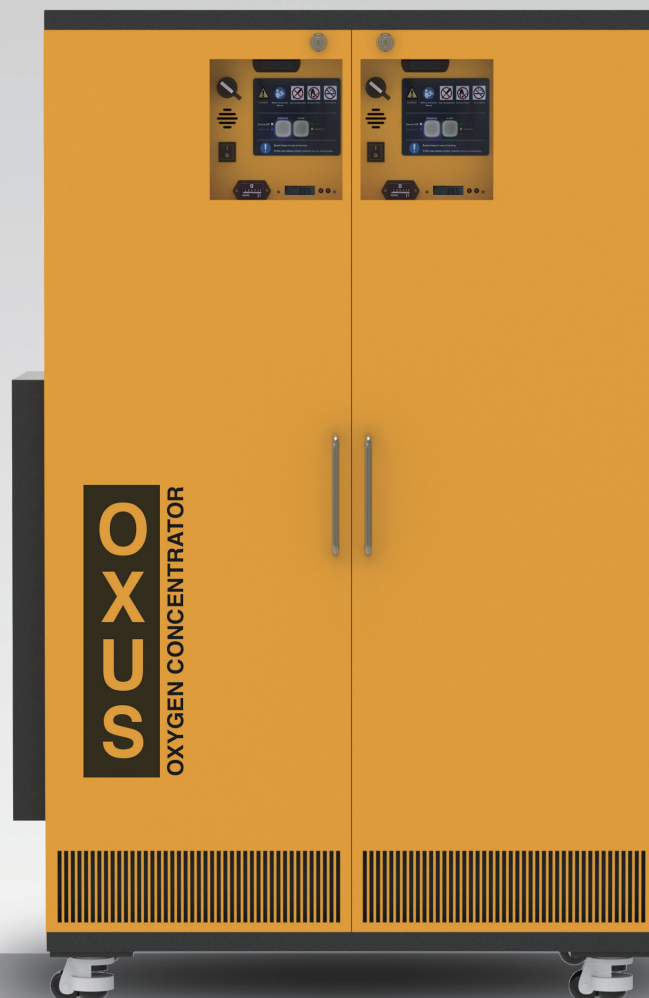

Product safety

- Power failure alarm function
- Low pressure alarm function

Oxygen supply efficiency

- Built-in booster(compressor)
- The product turns off automatically when the set pressure reaches 6.5 bar.


Maintenance differentiation

Even if one unit fails, the equipment does not turn off.

Medical oxygen supply system

As a system built to replace the continuous supply of liquefied and high-pressure oxygen to the central gas piping system used in medium, large hospitals, and nursing clinics medical oxygen supply system(RAK-Series) can continuously operate and supply oxygen to the central piping system for 24 hours, 365 days a year.

- Uninterrupted oxygen supply with continuous oxygen generation
- Excellent cost savings compared to normal liquid oxygen
- Space saving by eliminating the need for high-pressure vessel storage
- Saving maintenance costs by not using high-pressure cylinders

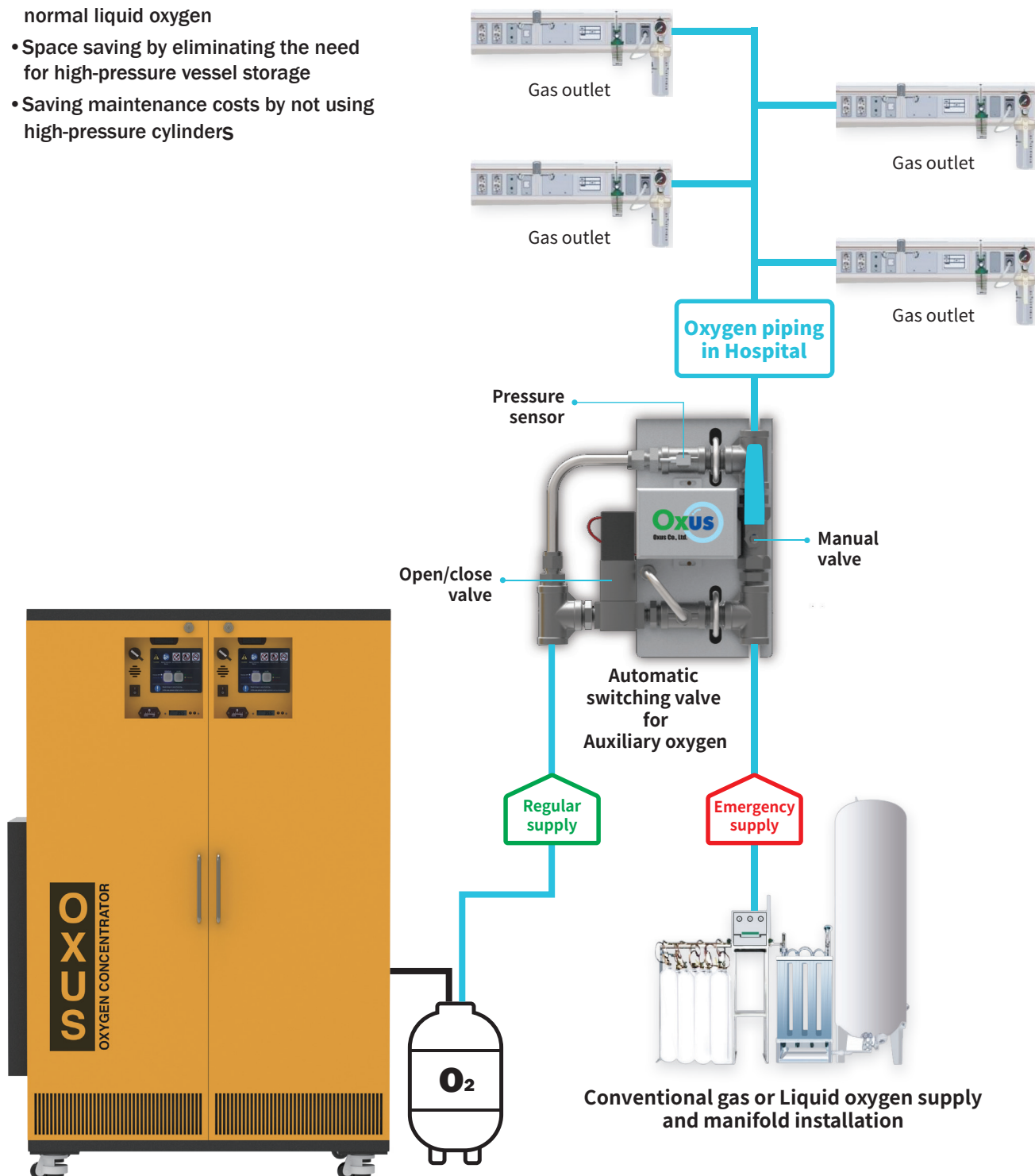




Why should you choose Oxus oxygen concentrator?


Pros of oxygen generator

- Ease of maintenance
- Cost effective
- Space saving

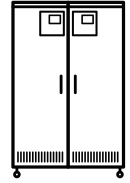

Problems with liquid and gas oxygen

- Hassle of oxygen tank filling and replacement
- Large space occupation
- Danger of high-pressure tank
- Danger of damaging facilities while moving tanks

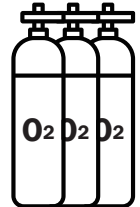

Product Detail

Since each units are connected in parallel, problems can be solved with simple unit replacement.




- 1,5 Control unit**
Operate the device and display the status
- 2,3,4,6,7,8 Oxygen unit**
Each unit produces oxygen



- 1 Control unit**
Operate the device and display the status
- 2,3,4 Oxygen unit**
Each unit produces oxygen

Units can be individually turned on/off depending on the amount of oxygen required.

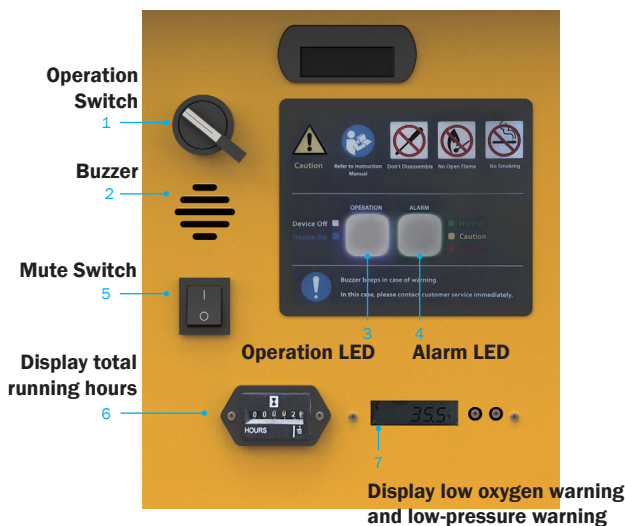

Optional Accessory

1. Oxygen Storage Tank

- This product requires Oxygen Storage Tank.
- Minimum tank volume is 100 Liter for 30LPM oxygen output.
- The picture illustrates 100 + 100 Liter tank arrangement (supplied by Oxus) for 60LPM Model.
- Any type of tank (oxygen rated above 8 bar) can be used for installation.

2. PC BOX

- Monitoring system possible





1 Day
↓

18
cylinders



RAK-U06M2E

Total production in one day: 86,400L

Gas oxygen (As for 40L, 120bar): 18 cylinders

Oxygen Flow	60±3 LPM
Oxygen Purity	93±3%
Oxygen Pressure	Max. 6.5 bar
Voltage(Frequency)	230V (50Hz)
Power Consumption	4.25 KW
Operating Environment	0°C ~ 40°C, 0% ~ 90 %, 1090hPa~800hPa
Noise Level	Max. 70dB
Weight	260kg
Dimension(W x D x H)	960 x 600 x 1,550mm
Electrical protection type	Class1



1 Day
↓

9
cylinders



RAK-U03M2E

Total production in one day: 43,200L

Gas oxygen (As for 40L, 120bar): 9 cylinders

Oxygen Flow	30±3 LPM
Oxygen Purity	93±3%
Oxygen Pressure	Max. 6.5 bar
Voltage(Frequency)	230V (50Hz)
Power Consumption	2.15 KW
Operating Environment	0°C ~ 40°C, 0% ~ 90 %, 1090hPa~800hPa
Noise Level	Max. 70dB
Weight	130kg
Dimension(W x D x H)	480 x 600 x 1,550mm
Electrical protection type	Class1